

# SAM Safety Seal Guedel

SAM

SCOTT AIRWAY MANAGEMENT

The **SAM Safety Seal** is designed to be used together with the **SAM Guedel** during bag mask ventilation with bearded, misshapen face and denture-less patients to assist in providing a better seal. The device is also used with any patient where aerosolisation of patient secretions during bag mask ventilation is a concern.

The device is easy to use. Simply invert the **SAM Safety Seal** wings into a gloved hand and place the guedel as per standard practice. Once the guedel is in situ, release the wings so that they lay onto the patient's face.

Once the bag mask ventilation process is completed, simply grab hold of the **SAM Safety Seal** wings and fold around the guedel into your gloved hand and discard into an appropriate receptacle.

Scan the QR Code to watch a video on how to use the **SAM Safety Seal** products.

## Product Specifications

Silicone - NON Latex

Body Height 2cm


SINGLE USE ONLY

Silicone Variant thickness - 1.5-2.5mm

Wing Dimensions 17cmx10cm

SAM SAFETY SEAL ARTG 331 999

SAM GUEDEL AIRWAY ARTG 353 530

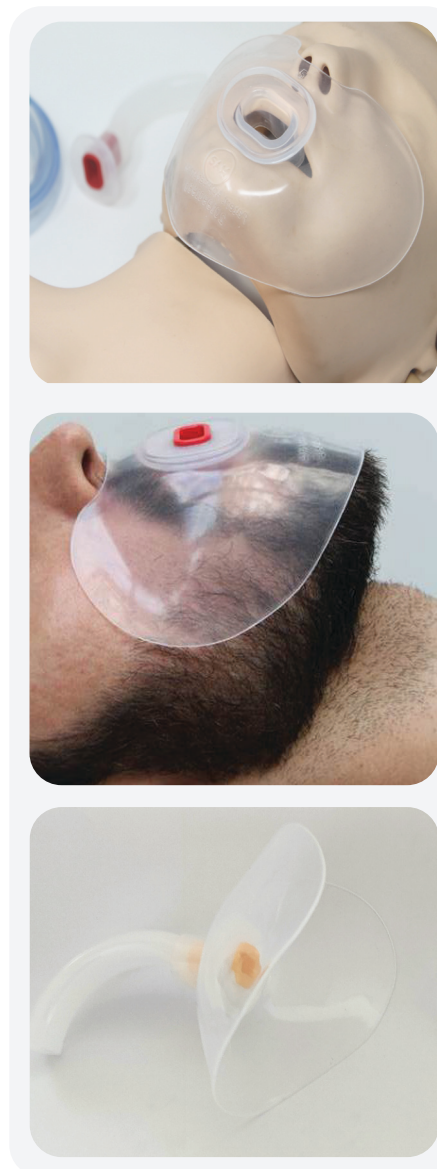

Product	Product Order Code	Quantity
SAM Safety Seal Guedel 9cm (Yellow)	SSSGN09	50
SAM Safety Seal Guedel 10cm (Red)	SSSGN10	50



**STATPACKS**  
PLUS 

**Raf Sciezka**

DISTRIBUTOR WA/SA

+61 400 990 465

raf.sciezka@statpacksplus.com.au