
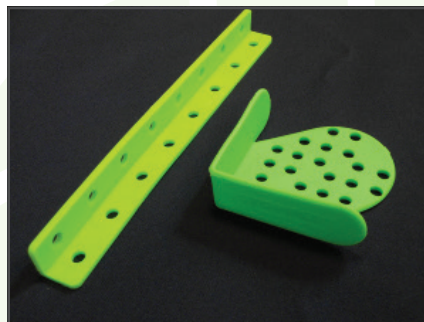


Problems with sterile wrap rips and tears? Try Stopper Care System Sterile Validated Corner Protectors.

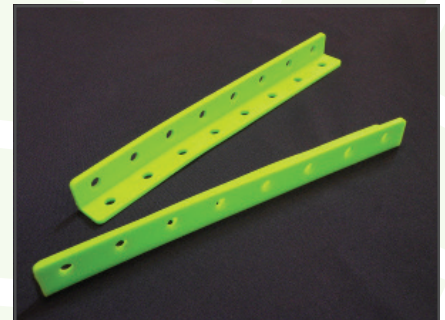

Validated:
Steam
ETO
Plasma

- Does not retain moisture
- Reduces amount of sterile wrap
- Will avoid tray pressing against wrap
- Saves OR and CSP time, money and frustration
- Reusable and economical

When the Stopper Care System is used on a surgical tray the surgical tray is completely elevated off the surgical wrapping.

The weight of the tray is totally transferred on to the Edge Protector and Stopper® Corner Shield, allowing for complete circulation of the sterilant and better air circulation for faster drying times.

In Steam Sterilization, this elevation also allows easier sliding of the trays since the weight of the tray is not sitting on top of the surgical wrappings. Since most of the instrumentation weight inside of the heavy surgical tray is concentrated in the middle, the tray tends to bow somewhat in the center of the tray.

This concentrated weight sits on the sterile wrapping and upon moving the tray (sliding) or pulling off a wire racked shelf (storage) the weight of the heavy tray could cause rips, holes and tears, especially in the underside of the wrapped surgical tray.

By using the Stopper Care System and the weight of the surgical tray being elevated and supported, the above is less likely to occur.

The Stopper Care System is reusable with both products being made of the same tough material, recyclable and environmentally friendly.

- Provides complete surgical tray protection from holes, rips, and tears.
- Elevates the complete tray away from the blue wrap for sliding of the tray and helps eliminate holes, rips, and tears on the bottom.
- Helps prevent holes, rips, and tears along the horizontal edges of heavy, sharp surgical trays
- Elevation allows for better drying time and sterilant dispersion
- Can be used with round corner trays
- Does not retain moisture
- Can be used in Steam, Eto, and Plasma (Sterrad) Sterilization
- Edge protectors can be cut to size to fit any size surgical tray
- Low cost, Multi-use, Reusable products



**LATERAL
MEDICAL**

Contact: (07) 3139 1207 or
sales@lateralmedical.com